

Keysight Scanning Tunneling Microscopy Scanners

Data Sheet



Overview

Scanning tunneling microscopy (STM) is a scanning probe microscopy (SPM) imaging technique that takes advantage of the extreme distance sensitivity of the tunneling current between two conducting electrodes. By measuring the tunnel-current variations as a probe is scanned over a sample's surface, STM is able to deliver the highest resolution SPM images.

Keysight STM scanners are designed to deliver outstanding results on a variety of conducting materials. These low-current and ultra-low current STM scanners provide stable imaging at pico-ampere and sub-pico-ampere currents to resolve individual atoms and molecules. Over the past ten years, Keysight STM scanners have delivered superior research results and have consistently outperformed other STM scanners in achieving the best atomic resolution within the shortest time.

A hermetically sealed, top-down configuration provides complete isolation of the scanning elements and electronics from the imaging environment. This design allows total environmental control, fluid friendly operation, and superior thermal stability. As a result, samples can be imaged at high temperatures up to 250°C over an extended period of time (tested up to 10 hours).


Keysight's versatile STM scanners are available in two ranges, atomic (1 μ m) and small (10 μ m). Both scanners are capable of atomic resolution and are designed to work with Keysight's AFM/SPM microscope systems. In addition, an STM nose cone option is available for Keysight's multipurpose scanners.

Features and Benefits

- Wide range of current sensitivity provides guaranteed atomic resolution imaging of conducting surfaces
- Designed for imaging in ambient, controlled gas, or in a fluid environment
- Top-down configuration protects electronics and piezo elements from damage caused by harsh imaging environments
- Full compatibility with Keysight's modular AFM/SPM microscopes offers simple upgrade path for extended capabilities
- Easy fluid exchange permits greater EC-STM versatility

Temperature Control

While imaging with Keysight STM scanners, a wide range of sample temperatures (from -30°C to 250°C) can be precisely controlled with up to $\pm 0.025^\circ\text{C}$ accuracy, in ambient or in liquid.

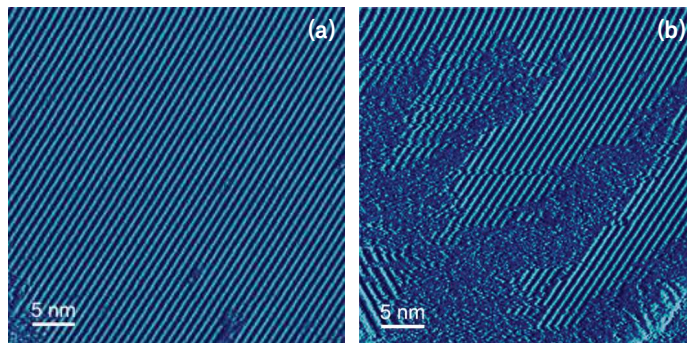

Scan size 40nm x 40nm

Figure 1. STM images of ordered bipyridine molecules undergo a phase transition with temperature change. (a) 29.7°C and (b) 33.51°C.

Application

Electrochemistry

Designed for electrochemistry, Keysight STM scanners deliver high resolution and ease of use in an extremely clean fluid environment.

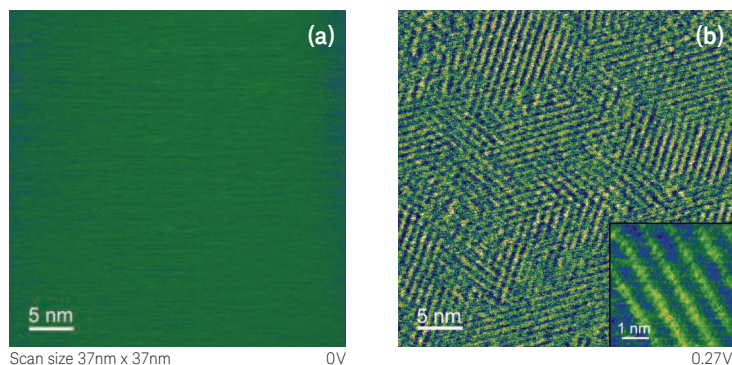


Figure 2. Potential-induced phase transition of 2,2'-bipyridine molecules on Au(111) surface imaged *in situ*. (a) At 0V (versus Ag/Ag+), randomly and loosely bound molecules were not observed. (b) At 0.27V, molecules were strongly bound to the surface and packed in ordered rows along the three distinct directions of the underneath [Au] atomic lattice. (b insert) High-resolution of closely packed individual bipyridine molecules at 0.27 V. 4.6nm x 4.6nm.

Molecular Chemical Identification

Keysight STM scanners have an integral role in the identity of individual molecules based on their chemical properties.

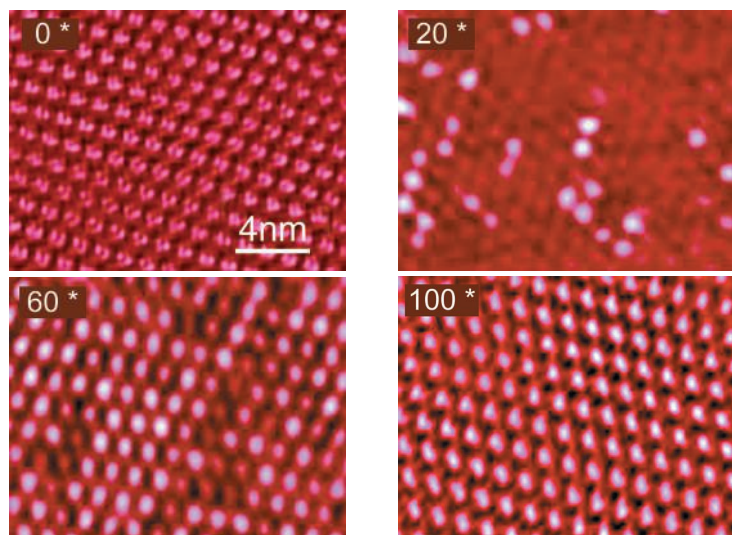


Figure 3. Protoporphyrin and Fe-Protoporphyrin on Au(111) imaged with EC-STM. 18nm x 14nm.

Specifications

Atomic STM

Scan range:	1 μm in XY, 0.7 μm in Z
Noise level:	<0.6 $\text{\AA}$ in XY, <0.06 $\text{\AA}$ in Z
Sensitivity:	0.1 nA/V, 1 nA/V (standard), or 10 nA/V
STM probe:	0.25mm Pt-Ir or W wire

Small STM

Scan range:	10 μm in XY, 1.6 μm in Z
Noise level:	<1 $\text{\AA}$ in XY, <0.6 $\text{\AA}$ in Z
Sensitivity:	0.1 nA/V, 1 nA/V (standard), or 10 nA/V
STM probe:	0.25mm Pt-Ir or W wire

AFM Instrumentation from Keysight Technologies

Keysight Technologies offers high precision, modular AFM solutions for research, industry, and education. Exceptional worldwide support is provided by experienced application scientists and technical service personnel. Keysight's leading-edge R&D laboratories are dedicated to the timely introduction and optimization of innovative, easy-to-use AFM technologies.

www.keysight.com/find/afm

For more information on Keysight Technologies' products, applications or services, please contact your local Keysight office. The complete list is available at: www.keysight.com/find/contactus

Americas

Canada	(877) 894 4414
Brazil	55 11 3351 7010
Mexico	001 800 254 2440
United States	(800) 829 4444

Asia Pacific

Australia	1 800 629 485
China	800 810 0189
Hong Kong	800 938 693
India	1 800 112 929
Japan	0120 (421) 345
Korea	080 769 0800
Malaysia	1 800 888 848
Singapore	1 800 375 8100
Taiwan	0800 047 866
Other AP Countries	(65) 6375 8100

Europe & Middle East

Austria	0800 001122
Belgium	0800 58580
Finland	0800 523252
France	0805 980333
Germany	0800 6270999
Ireland	1800 832700
Israel	1 809 343051
Italy	800 599100
Luxembourg	+32 800 58580
Netherlands	0800 0233200
Russia	8800 5009286
Spain	0800 000154
Sweden	0200 882255
Switzerland	0800 805353

	Opt. 1 (DE)
	Opt. 2 (FR)
	Opt. 3 (IT)
United Kingdom	0800 0260637

For other unlisted countries:
www.keysight.com/find/contactus

This information is subject to change without notice.
© Keysight Technologies, 2009–2014
Published in USA, August 26, 2014
5989-5997EN
www.keysight.com